

# QBot Platform

Play

© 2024 Quanser Inc., All rights reserved.

Quanser Inc.  
119 Spy Court  
Markham, Ontario  
L3R 5H6  
Canada

info@quanser.com  
Phone: 19059403575  
Fax: 19059403576  
Printed in Markham, Ontario.



For more information on the solutions Quanser Inc. offers, please visit the web site at:  
<http://www.quanser.com>

This document and the software described in it are provided subject to a license agreement. Neither the software nor this document may be used or copied except as specified under the terms of that license agreement. Quanser Inc. grants the following rights: a) The right to reproduce the work, to incorporate the work into one or more collections, and to reproduce the work as incorporated in the collections, b) to create and reproduce adaptations provided reasonable steps are taken to clearly identify the changes that were made to the original work, c) to distribute and publicly perform the work including as incorporated in collections, and d) to distribute and publicly perform adaptations. The above rights may be exercised in all media and formats whether now known or hereafter devised. These rights are granted subject to and limited by the following restrictions: a) You may not exercise any of the rights granted to You in above in any manner that is primarily intended for or directed toward commercial advantage or private monetary compensation, and b) You must keep intact all copyright notices for the Work and provide the name Quanser Inc. for attribution. These restrictions may not be waved without express prior written permission of Quanser Inc.

# QBot Platform – Application Guide

## Play

### Why play?

---

Before getting into the lab content using the QBot Platform, it is important to first get a good understanding of how to interface with the robot, deploy code and models to it, and be able to understand the available hardware capabilities. This first lab is a chance for you to play around with the sensors, camera, lidar, and driving capabilities of the robot to get a good feel for how it works!

## QBot Platform Overview

The QBot Platform is equipped with a variety of different hardware: two motors, corresponding encoders (position) and tachometers (speed) for each motor, a downward facing grayscale global shutter camera, a forward-facing colour and depth camera, a lidar, user programmable LEDs and an LCD display. The QBot Platform also includes magnetic attachment points on the top surface and attachment slots for adding your own hardware.

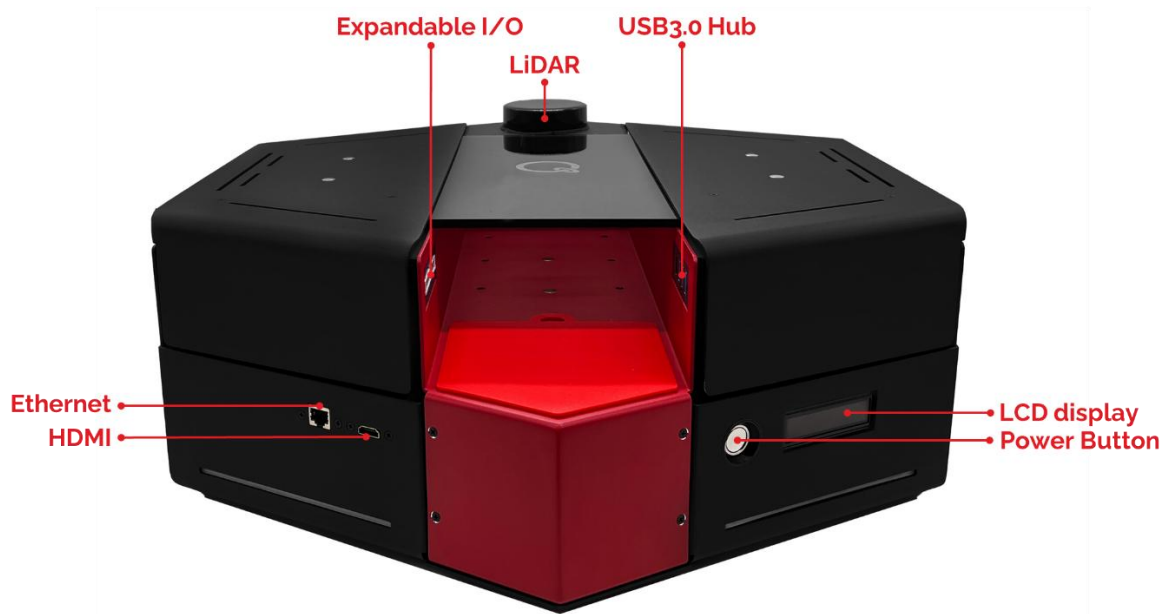


Figure 1: QBot Platform Back View

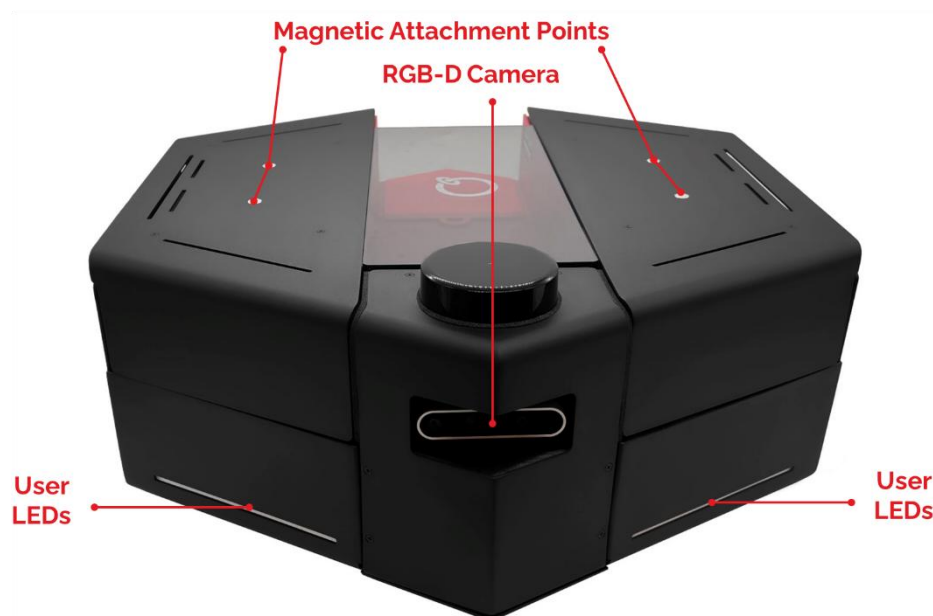


Figure 2: QBot Platform Front View

## Getting Started

The QBot Platform is a differential drive robot designed to allow students to understand key fundamental topics in robotics: sensors and actuators, control systems, programming and software, robot perception, path planning and mapping and the list goes on!

Before you begin this lab, make sure that the following criteria have been met:

- The QBot Platform has been setup and tested.
  - o See the **QBot Platform Quick Start Guide** for details on this step.
  - o Review the **QBot Platform User Manuals** for more information.
- You are familiar with Python or Simulink.

When you are ready, start the [Lab 1 Play – Lab Procedure](#).