

Lab Procedure

Pendulum Interfacing & Inertia Virtual

Introduction

Ensure the following:

1. You have reviewed the [Application Guide – Pendulum Interfacing and Inertia](#)
2. Make sure you have Quanser Interactive Labs open in the Qube 3 – Pendulum → Pendulum Workspace.
3. Launch MATLAB and browse to the working directory that includes the Simulink models for this lab.

The **Hardware Interfacing** and **Filtering** labs explained the basic blocks to read and write from the Qube-Servo 3. For simplicity, all labs forward will use a Qube-Servo 3 block that sets up the system beforehand and outputs the available information from the Qube.

This lab will focus on understanding the pendulum attachment and will be separated in three sections. The first one will help you develop the pendulum conventions of movement. The last two will be focused on calculating and measuring the moment of inertia of the pendulum.

Developing the Pendulum Conventions

Using the gains found to convert tachometer counts/s into rads/s from the instrumentation labs, update the gains in the [qs3_pendulum_interfacing.slx](#) file. This model should read the base motor shaft and pendulum encoder positions of the Qube-Servo 3 and apply a small voltage to the motor. The model will then be modified to adhere to the conventions of the pendulum model, as laid out in the concept review.

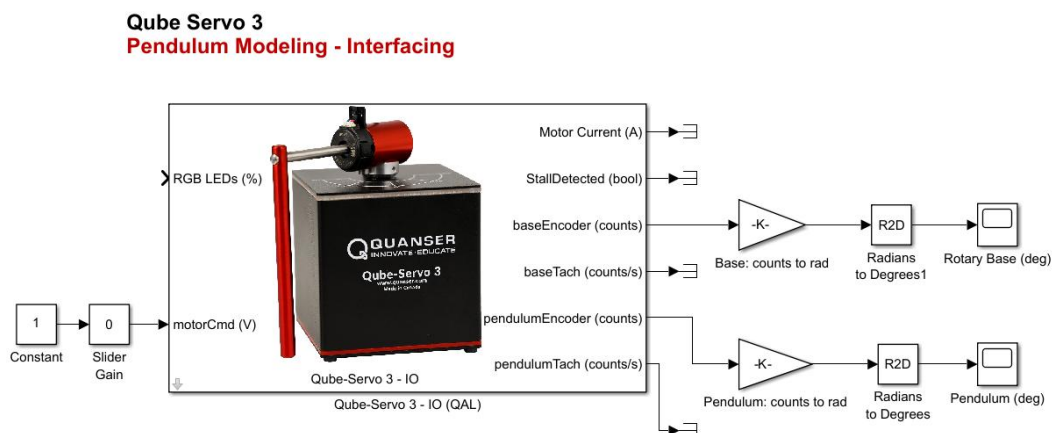



Figure 1: Simulink model used to interface with the Qube-Servo 3 using the pendulum module.

Setup the model as follows:

1. Setup the encoder gains on **baseEncoder** and **pendulumEncoder** to read the angles in radians.
2. Connect the measured angles to scopes but display them in degrees. You can do this adding the **Radians to Degrees** block into your model, by double clicking on an empty part of your model and typing the block name or through another **Gain** block.
3. Ensure the voltage sent to the Qube-Servo 3 can be adjusted. This can be done through a **Constant** block and a **Slider Gain** block. Make sure the voltage sent to the motor is limited to $\pm 1V$.
4. Run the QUARC controller using the Run  button on the **Simulation** tab.
5. Click Lift Pendulum in the virtual workspace. Examine the response in the scopes based on the side the pendulum lifts and the way the base moves once it starts falling.



6. If the rotations performed in the previous step do not match the positive rotation conventions of the pendulum model, modify the Simulink model so that they do.
Note: The Simulink model may need to be stopped and have **Gain** blocks added so that CCW is the positive direction.
7. Apply a small voltage (0.5V) to the motor by using the **Slider Gain** block or modifying the **Constant** block.
8. If the voltage applied does not adhere to the expected rotation from a positive voltage, modify the Simulink model so that it does.
9. Stop the Simulink model.
10. Add a **Modulus** and **Bias** block to the Simulink model so that the pendulum angle α is 0 rad in the upright position and is bounded by $\pm\pi$ rad. The model should look similar to Figure 2 below:

Qube Servo 3 Pendulum Modeling - Interfacing

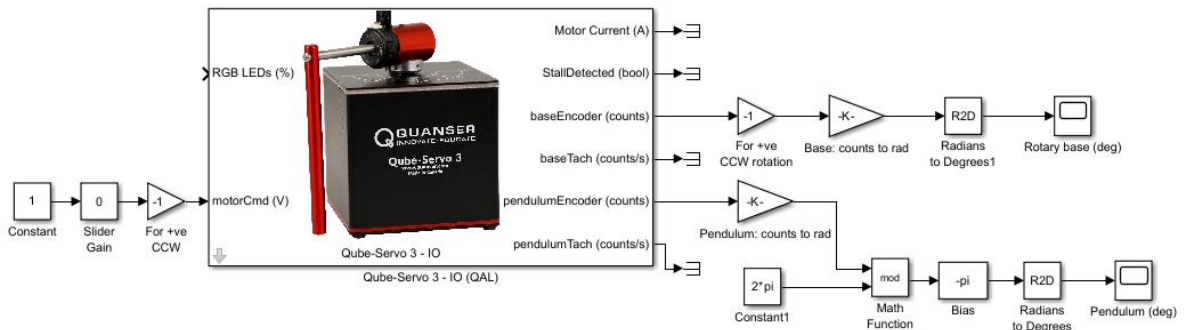


Figure 2: Simulink model that includes the modulus and bias blocks.

11. Run the QUARC controller using the Run button on the Simulation tab.
12. Rotate the pendulum to the upright position to ensure the angle is 0 rad using the lift pendulum button. Capture a screenshot of the Pendulum (deg) scope demonstrating the pendulum moving from the downright position to the upright position.
13. Once all conventions have been validated, stop the Simulink model.
14. Take a screenshot of the completed Simulink model with all the changes made.
15. Close your model. Ensure you save a copy of the files for review later.

Calculating the Moment of Inertia

The goal is to derive the moment of inertia of the pendulum module on the Qube-Servo 3 experimentally and analytically. The model [qs3_moment_of_inertia.slx](#) measures the pendulum angle using the encoder and converts from encoder counts to degrees.

Qube Servo 3 Pendulum Modeling - Moment of Inertia

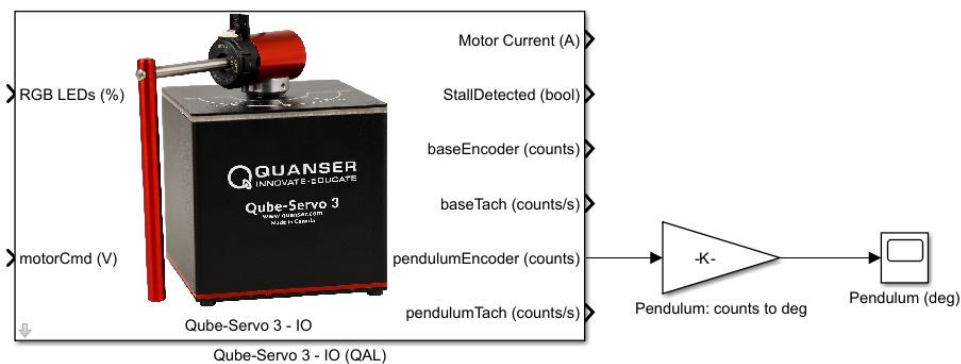


Figure 3: Simulink model that displays the measured pendulum angle.

16. Find the moment of inertia acting about the pendulum using the free-body diagram from the concept review. Use the parameters as defined in the Qube-Servo 3 User Manual. Use the file `qs3_rotpen_param.m` to fill out your values and do the calculations.

Hint: For solid objects with uniform density, you can express the differential mass in terms of differential length.

Measuring the Moment of Inertia

17. Run the QUARC controller using the Run button on the Simulation tab.
18. On the Settings in Quanser Interactive labs, set the Lock Servo base to on and then click lift pendulum. This will make sure the pendulum swings freely without the base rotating.
19. Measure the number of cycles (n_{cyc}) and the duration of the cycles (Δt) from the pendulum position scope. Your response should look similar to Figure 4:

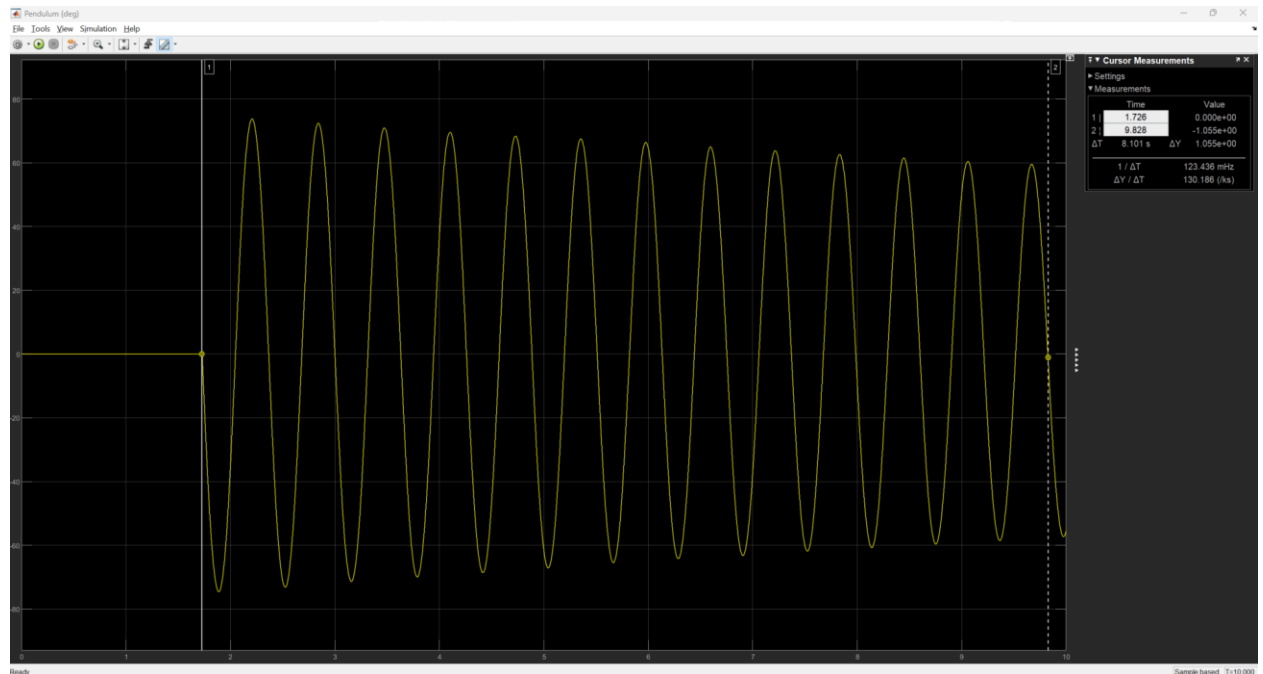


Figure 4: Sample response of perturbed pendulum.

20. Calculate the moment of inertia for the pendulum using n_{cyc} and Δt .
21. Once the calculations have been validated, close your model. ensure you save a copy of the files for review later.
22. Close Quanser Interactive Labs.