

# Lab Procedure

## Parameter Estimation

### Introduction

Ensure the following:

1. You have reviewed [Application Guide - Parameter Estimation](#)
2. The Qube-Servo 3 has been previously tested, is ON and connected to the PC.
3. Inertia disc load is attached to the Qube-Servo 3.
4. Launch MATLAB and browse to the working directory that includes the Simulink models for this lab.

The **Hardware Interfacing** and **Filtering** labs explained the basic blocks to read and write from the Qube-Servo 3. For simplicity, all labs forward will use a Qube-Servo 3 block that sets up the system beforehand and outputs the available information from the Qube.

Using the gains found to convert tachometer counts/s into rads/s from the instrumentation labs, use the [qs3\\_parameter\\_estimation.slx](#) file to design a model that applies a constant voltage to the motor and reads the motor current of the Qube-Servo 3 as shown in Figure 1.

### Qube Servo 3 Modeling - Parameter Estimation

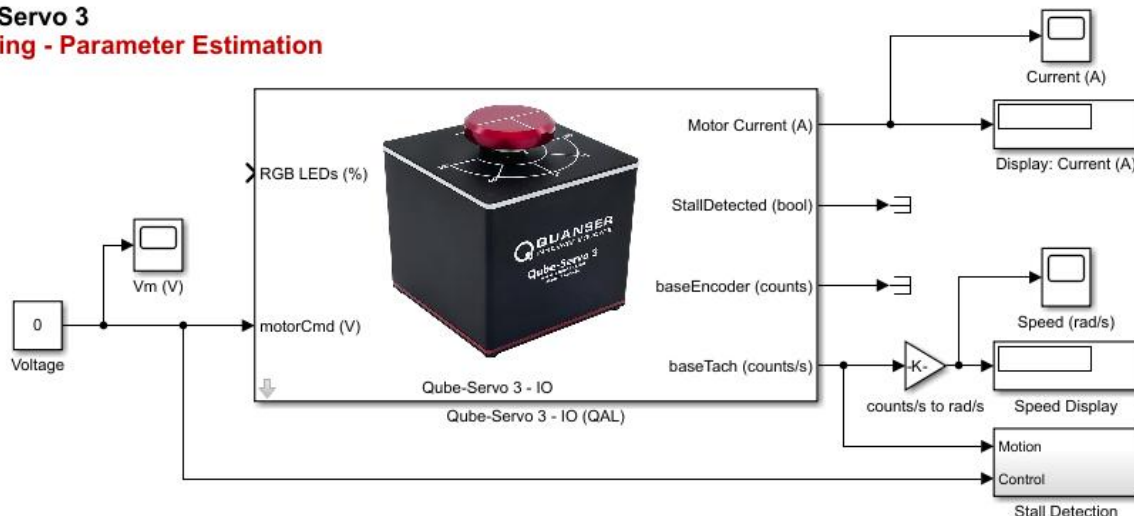


Figure 1: Model applying a voltage to the Qube and measuring motor current and speed.

During the lab, you will be experimentally estimating the motor resistance  $R_m$ . This can be done by applying constant voltages to the motor and measuring the corresponding current *while holding the motor shaft stationary*. You will also be experimentally estimating the motor torque constant  $k_m$ . Refer to [Concept Review – Motor Equations](#).

## Motor Resistance Estimation

1. Derive an expression that will allow you to solve for the motor resistance  $R_m$  by applying constant voltages to the motor and stalling it and measuring the corresponding current.
2. Ensure the inertial disc load is mounted on the Qube-Servo 3.
3. Build and run the QUARC controller using the **Monitor & Tune**  button on the **Hardware** or **QUARC** tab.
4. To experimentally estimate the motor resistance, apply a set of voltages to the Qube-Servo 3 using the model in Figure 1. For each measurement, *hold the motor shaft stationary* by grasping the inertial disc load to stall the motor. Record the current measurement displayed in the Current (A) display. To prevent any issues with the Qube-Servo, please do not hold the motor stalled for more than 2 seconds at a time.
5. If you accidentally hold the motor stalled for longer at higher voltages, the servo will stop spinning. If that happens, stop the model and start it again.
6. Fill the following table with the measured current for different voltage and calculate the corresponding resistance from the equation you derived in step 1.

| Applied Voltage<br>$V_m$ (V) | Measured Current<br>$I_m$ (A) | Resistance<br>$R_m$ ( $\Omega$ ) |
|------------------------------|-------------------------------|----------------------------------|
| -5                           |                               |                                  |
| -4                           |                               |                                  |
| -3                           |                               |                                  |
| -2                           |                               |                                  |
| -1                           |                               |                                  |
| +1                           |                               |                                  |
| +2                           |                               |                                  |
| +3                           |                               |                                  |
| +4                           |                               |                                  |
| +5                           |                               |                                  |

Table 1. Motor Resistance Experimental Results

7. Take the average of all the measure resistance values and compare this with the motor resistance value from the Qube-Servo 3 User Manual.
8. Stop your model.

## Motor Back-EMF Estimation

9. Derive an expression that will allow you to solve for the motor torque constant  $k_m$  by applying constant voltages to the motor and measuring the steady state current and speed (in rad/s). Assume you already know the motor resistance  $R_m$ .
10. Build and run the QUARC controller using the **Monitor & Tune**  button on the **Hardware** or **QUARC** tab.
11. To experimentally estimate the motor back-EMF constant, repeat the same procedure by applying different voltage to the Qube-Servo 3 with the motor free to spin (i.e. do not stall the motor) and record the measured speed and current in Table 2.

| Applied Voltage<br>$V_m$ (V) | Measured Speed<br>$\omega_m$ (rad/s) | Measured Current<br>$I_m$ (A) | Back-EMF<br>$k_m$ (V-s/rad) |
|------------------------------|--------------------------------------|-------------------------------|-----------------------------|
| -5                           |                                      |                               |                             |
| -4                           |                                      |                               |                             |
| -3                           |                                      |                               |                             |
| -2                           |                                      |                               |                             |
| -1                           |                                      |                               |                             |
| +1                           |                                      |                               |                             |
| +2                           |                                      |                               |                             |
| +3                           |                                      |                               |                             |
| +4                           |                                      |                               |                             |
| +5                           |                                      |                               |                             |

Table 2. Motor Resistance Experimental Results

12. Take the average of the measured back-EMF values and compare this with the motor back-EMF value from the Qube-Servo 3 User Manual.
13. The motor shaft of the Qube-Servo 3 is attached to a *load hub* and a *disk load*. Based on the parameters given in the Qube-Servo 3 User Manual, calculate the equivalent moment of inertia that is acting on the motor shaft. Keep note of your procedure as it may be required in the assessment.
14. Formulate the differential equation for  $\omega_m$  using the equations from the concept review outlined in the application guide. (*Hint:* Obtain the Voltage  $V_m(s)$  to Speed  $\Omega_m(s)$  transfer function by applying a Laplace Transform to the derived differential equation.) Keep note of your procedure as it may be required in the assessment.
15. Stop and close your model.
16. Power OFF the Qube-Servo 3.