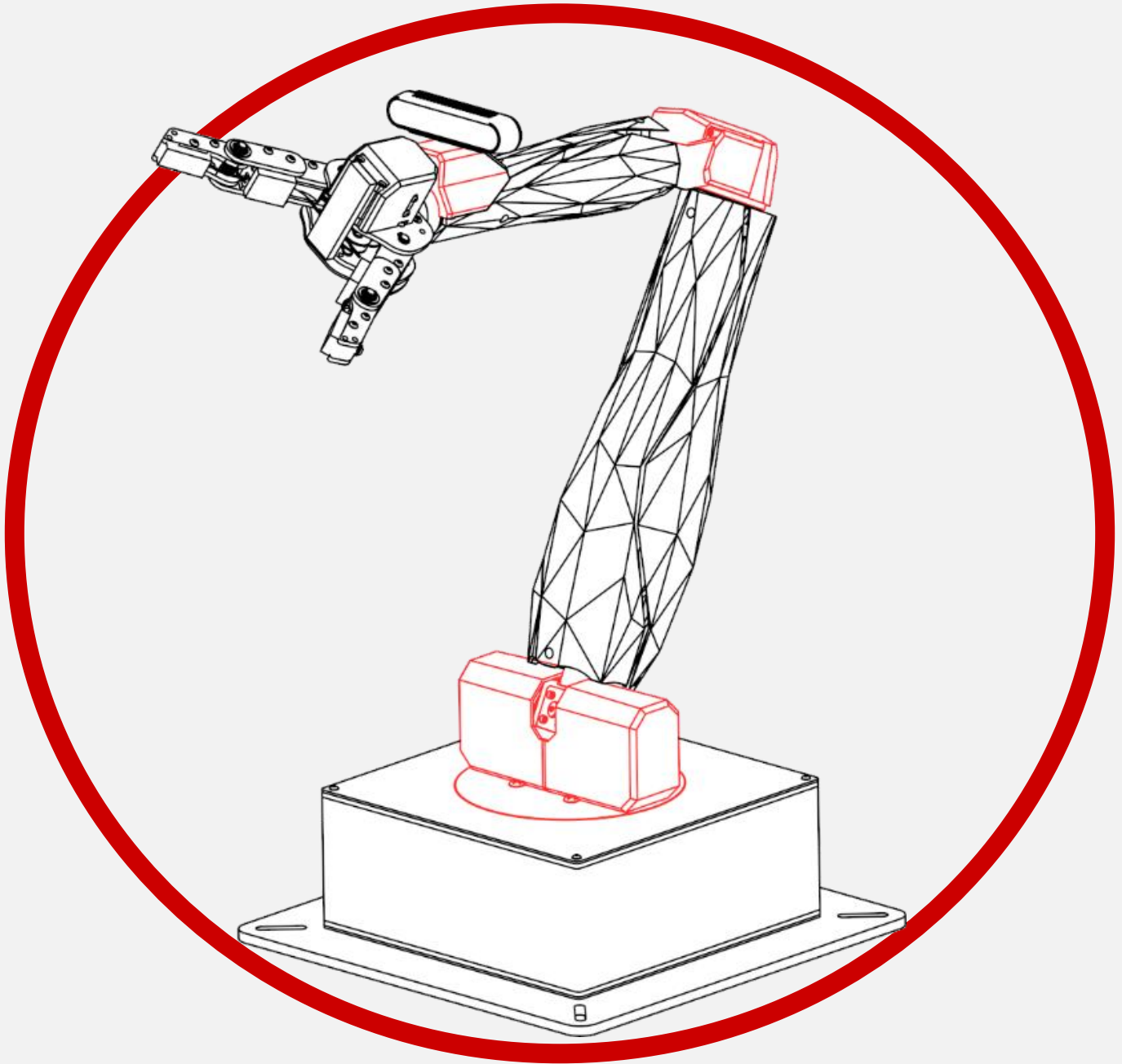




QArm

Tool Manipulation

V1.0 – 13th February 2025

© 2025 Quanser Consulting Inc., All rights reserved.
For more information on the solutions Quanser offers,
please visit the web site at: <http://www.quanser.com>



Quanser Consulting Inc. info@quanser.com
119 Spy Court Phone : 19059403575
Markham, Ontario Fax : 19059403576
L3R 5H6, Canada printed in Markham, Ontario.

This document and the software described in it are provided subject to a license agreement. Neither the software nor this document may be used or copied except as specified under the terms of that license agreement. Quanser Consulting Inc. ("Quanser") grants the following rights: a) The right to reproduce the work, to incorporate the work into one or more collections, and to reproduce the work as incorporated in the collections, b) to create and reproduce adaptations provided reasonable steps are taken to clearly identify the changes that were made to the original work, c) to distribute and publicly perform the work including as incorporated in collections, and d) to distribute and publicly perform adaptations. The above rights may be exercised in all media and formats whether now known or hereafter devised. These rights are granted subject to and limited by the following restrictions: a) You may not exercise any of the rights granted to You in above in any manner that is primarily intended for or directed toward commercial advantage or private monetary compensation, and b) You must keep intact all copyright notices for the Work and provide the name Quanser for attribution. These restrictions may not be waived without express prior written permission of Quanser.



Caution

This equipment is designed to be used for educational and research purposes and is not intended for use by the public. The user is responsible for ensuring that the equipment will be used by technically qualified personnel only. Users are responsible for certifying any modifications or additions they make to the default configuration.

QArm – Application Guide

Tool Manipulation

Why explore Tool Manipulation?

In the previous labs, we analyzed the forward and inverse kinematic formulations of a serial manipulator, which related the pose of the end-effector to its joint variables. The speed of the manipulator's joints or end-effector was not considered. For a task such as spray painting, not only do we care about precision in terms of 'where' the manipulator is painting, but also the uniformity with which it paints. To regulate a uniform speed for the end-effector, we must consider the relationship between the joint velocities of the manipulator and the task space velocities of the end-effector. This relationship, called the differential kinematics formulation, involves taking the time derivative of the forward kinematics formulation.

For the previous labs, the low-level joint controllers moved the manipulator to the desired joint-space position. Here, the controller regulates the arm based on a joint space error but does not monitor the task space performance. To regulate task space motion, a state machine is introduced that automates decision making. The manipulator will be commanded to move towards a desired setpoint at a constant rate until it gets close enough, regulating both position and speed.

Differential Kinematics

The forward kinematics formulation that related the joint angles to the end-effector position and orientation are outline below. (Recall that c_1 stands for $\cos \theta_1$, s_1 stands for $\sin \theta_2$, s_{23} stands for $\sin(\theta_2 + \theta_3)$ etc.).

$$\begin{aligned} P &= f_{FPK}(\Theta) \\ p_x &= \lambda_2 c_1 c_2 - \lambda_3 c_1 s_{23} \\ p_y &= \lambda_2 s_1 c_2 - \lambda_3 s_1 s_{23} \\ p_z &= \lambda_1 - \lambda_2 s_2 - \lambda_3 c_{23} \\ \gamma &= \theta_4 \end{aligned} \quad (1)$$

where Θ is the vector of the joint states,

$$\Theta = \begin{bmatrix} \theta_1 \\ \theta_2 \\ \theta_3 \\ \theta_4 \end{bmatrix} \quad (2)$$

and p is the position of the end-effector [4] (expressed in base frame {0})

$$p = \begin{bmatrix} p_x \\ p_y \\ p_z \end{bmatrix} = {}^0p_4 \quad (3)$$

To formulate the differential kinematics, we differentiate the equation set 1 with respect to time

$$\begin{aligned} v_x &= \frac{dp_x}{dt} = -\lambda_2 s_1 c_2 \dot{\theta}_1 - \lambda_2 c_1 s_2 \dot{\theta}_2 + \lambda_3 s_1 s_{23} \dot{\theta}_1 - \lambda_3 c_1 c_{23} \dot{\theta}_2 - \lambda_3 c_1 c_{23} \dot{\theta}_3 \\ v_y &= \frac{dp_y}{dt} = \lambda_2 c_1 c_2 \dot{\theta}_1 - \lambda_2 s_1 s_2 \dot{\theta}_2 - \lambda_3 c_1 s_{23} \dot{\theta}_1 - \lambda_3 s_1 c_{23} \dot{\theta}_2 - \lambda_3 s_1 c_{23} \dot{\theta}_3 \\ v_z &= \frac{dp_z}{dt} = -\lambda_2 c_2 \dot{\theta}_2 + \lambda_3 s_{23} \dot{\theta}_2 + \lambda_3 s_{23} \dot{\theta}_3 \\ \dot{\gamma} &= \frac{d\gamma}{dt} = \dot{\theta}_4 \end{aligned} \quad (4)$$

Grouping the terms based on joint velocities and representing it in matrix form yields,

$$\begin{aligned} V &= {}^0J \dot{\Theta} \\ \begin{bmatrix} v_x \\ v_y \\ v_z \\ \dot{\gamma} \end{bmatrix} &= \begin{bmatrix} -\lambda_2 s_1 c_2 + \lambda_3 s_1 s_{23} & -\lambda_2 c_1 s_2 - \lambda_3 c_1 c_{23} & -\lambda_3 c_1 c_{23} & 0 \\ \lambda_2 c_1 c_2 - \lambda_3 c_1 s_{23} & -\lambda_2 s_1 s_2 - \lambda_3 s_1 c_{23} & -\lambda_3 s_1 c_{23} & 0 \\ 0 & -\lambda_2 c_2 + \lambda_3 s_{23} & \lambda_3 s_{23} & 0 \\ 0 & 0 & 0 & 1 \end{bmatrix} \begin{bmatrix} \dot{\theta}_1 \\ \dot{\theta}_2 \\ \dot{\theta}_3 \\ \dot{\theta}_4 \end{bmatrix} \end{aligned} \quad (5)$$

where 0J is the manipulator Jacobian expressed in base frame {0}.

Background

The QArm content contains 3 labs that focus on velocity manipulation. They focus on tool manipulation, singularity identification and singularity avoidance respectively. This lab focuses on understanding differential kinematics to get task space velocities of the end effector.

Getting started

The goal of this lab is to understand how to avoid singularity positions based on the singularity identification lab before.

Ensure you have completed the following labs

- **Kinematic Manipulation Labs**

Before you begin this lab, ensure that the following criteria are met.

- The QArm has been setup and tested. See the QArm Quick Start Guide for details on this step.
- You are familiar with the basics of Simulink. See the [Simulink Onramp](#) for more help with getting started with Simulink.